



Deskset 17 & Deskset 19

ONYX Deskset Models at a Glance:

- Self-contained and easy to carry
- Magnification up to 85x (Deskset 17) and 95x (Deskset 19)
- Distance, document, and self views
- Camera mounts on left or right side
- 350 degree camera rotation
- Adjustable monitor height and tilt

Power and Productivity

ONYX Deskset models are ideal "on-the-go" low vision solutions for clear, powerful magnification that is portable, adaptable, and easy to use. They fit in any classroom, business, or home environment and move from one viewing situation to the next as quickly as you do.

ONYX Deskset models feature the versatile ONYX camera mounted on a rotating arm to a lightweight flat panel monitor. They provide you with the ultimate in distance and up-close magnification in a compact, easy-to-carry, self-contained unit that plugs into any standard wall socket.




Clear Images – Near and Far

ONYX Deskset models produce precise magnification in three distinct viewing modes. In the distance mode, even objects from across large rooms or auditoriums can be seen in sharp focus. The document mode magnifies items such as reading materials, product labels, and craft projects. With the self view, you get a magnified true mirror image close up. The camera remembers its settings in each position, so you can switch back and forth between viewing modes without refocusing or readjusting the camera.

Versatility for any Situation

The camera's mounting arm has three axes, and the camera itself rotates up to 350 degrees, enabling you to place it in an almost endless number of positions.

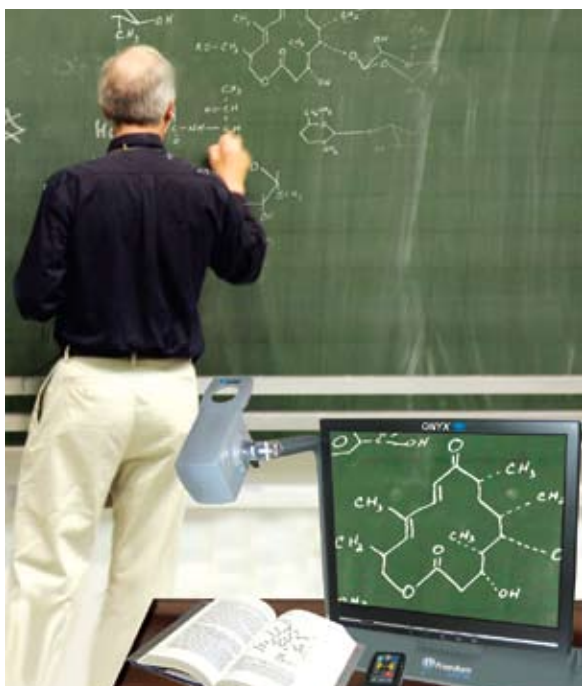

ONYX™

Freedom Scientific

Deskset 17 & Deskset 19

The camera can be mounted on either side of the monitor, and the monitor can be raised and lowered and tilted forward and backward for maximum comfort.

The many key features of ONYX Deskset models include:

- A choice of full color, enhanced black on white positive, enhanced white on black negative, or an inverse color mode that can be used to enhance contrast without loss of fine details
- A full range of brightness adjustments to eliminate glare and enable viewing in low light
- Masking and Reading Lines

- A compact, lightweight design with an integrated handle. The camera folds away behind the monitor, and the entire unit slips into the included rolling case
- The Find Feature easily establishes viewing perspective on an image prior to zooming in
- Focus Lock keeps an object or document crisp and clear when working or writing under the camera
- Freeze Frame provides ample time to carefully review images
- A unique, pocket-sized RF remote control
- A universal 100-240 volt power supply



Easy to take with you. Just grasp the convenient carrying handle, and slip your ONYX Deskset into its case. Expand the telescoping handle, and roll the case behind you.

Specifications

Weight:

- ONYX Deskset 17: 12 lbs / 5.44 kg
- ONYX Deskset 19: 15.4 lbs / 7 kg

Magnification at 11 3/8 in / 289 mm:

- ONYX Deskset 17: 85x
- ONYX Deskset 19: 95x

Camera-only weight: 14.8 oz / 420 g

Magnification at 3 ft / 91.44 cm

- ONYX Deskset 17: 21.3x
- ONYX Deskset 19: 23.8x

Working space under camera: 11 3/8 in / 289 mm

About Freedom Scientific

Freedom Scientific is the world's leading manufacturer of assistive technology products for those who are blind or vision impaired. For more information, e-mail info@freedomscientific.com

Designed & Assembled
in the USA 