



Swing-arm PC Edition



ONYX Swing-arm PC Edition at a Glance:

- Turn any PC into a powerful vision and image management tool
- Switch easily between magnified PC view and camera view
- Capture multiple magnified images
- Magnify up to 75x
- MAGic screen magnification software with speech included!

Power and Productivity

Bring your digital world and physical world together seamlessly in an intuitive work environment. This powerful new vision tool lets you work with a magnified PC screen and view magnified camera images - at the same time!


Featuring a small but powerful camera that provides magnification up to 75 times on a 15-inch

PC screen, the ONYX Swing-arm PC Edition creates clear magnified views of your documents and distant objects.

It's ideal for use in the workplace and classroom. The ONYX camera moves quickly and easily back and forth between three distinct views: distance, document, and mirror-image self views.

Complete View Management on Any PC

The ONYX Swing-arm PC Edition features USB connectivity to any PC and comes complete with Freedom Scientific's MAGic[®] screen magnification software with speech. The addition of new MAGic EyeMerge[™] features means you get powerful distance and up close magnification,

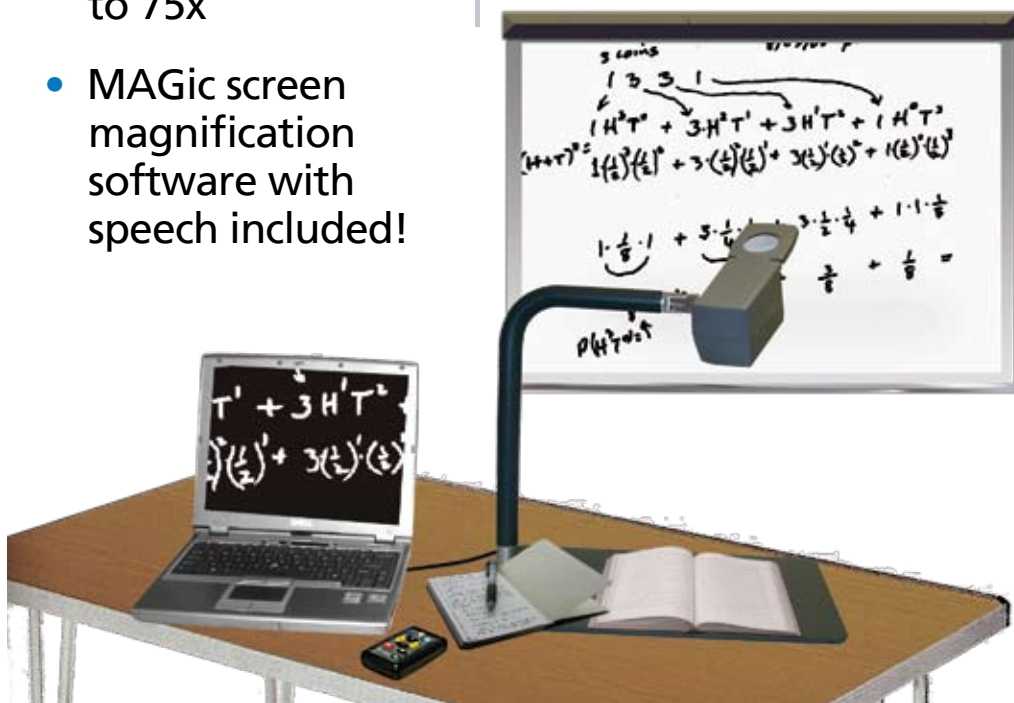

PC screen magnification, and the ability to customize, save and manipulate images to maximize productivity when working with a PC, documents, and/or viewing a presentation or lecture.

Productivity in Overdrive

Unique software features let you capture multiple snapshot views of a blackboard, for example, so that you can easily flip through them on your PC, rather than constantly panning and readjusting the camera. Capture blackboard or presentation images before they are erased or removed, and save them as your notes.

The many features of the ONYX Swing-arm PC Edition include:

- Find Feature
- Focus Lock
- Freeze Frame
- RF remote control
- Detachable camera
- Universal 110-240 volt power supply

- A compact, lightweight design that sets up and disassembles quickly for easy portability
- Camera can be powered from a PC via two USB ports
- MAGic screen magnification software with speech and MAGic EyeMerge™ features:
 - Switch easily between magnified PC view and camera view
 - Create multiple split screen views
 - Change text and background colors
 - Capture images
 - Flip easily between multiple snapshots
 - Adjust contrast levels

Specifications

Weight:	Magnification at
6 lbs / 2.72 kg	11 3/8 in / 289 mm:
Camera-only weight:	75x *
14.8 oz / 6.71 kg	Magnification at
Working space under	3 ft / 99.44 cm:
camera:	18.75x*
11 3/8 in / 289 mm	

Required: Support for Microsoft® DirectX® 9

*Based on 15" PC display

Designed & Assembled
in the USA 

About Freedom Scientific

Freedom Scientific is the world's leading manufacturer of assistive technology products for those who are blind or vision impaired. For more information, e-mail info@freedomscientific.com